FIND WELL-BEING

# USER MANUAL

## CONTACT US

Aukey International Limited

**Web:** [www.naipocare.com](http://www.naipocare.com)

**E-mail:** [wholesale@naipocare.com](mailto:wholesale@naipocare.com)

**Add.:** Room 601 6/F Tung Hip Commercial Building,  
248 Des Voeux Road Central Sheung Wan, Hong Kong

**HQ:** No. 102, Building P09, China South City, No. 1 Huanan Avenue,  
Pinghu Street, Longgang District, Shenzhen, China

## MADE IN CHINA

[/naipocare](#) [/naipocare](#) [/naipocare](#)

## oCuZen S

NECK & FULL BACK MASSAGER

# CONTENTS |

English

02

---

Deutsch

12

---

Français

22

---

Español

32

---

Italiano

42

---

## IMPORTANT SAFETY INSTRUCTIONS

When using the product, all the following instructions and basic precautions should be taken.

Please keep this manual for future reference over the duration of ownership.

## DANGER

To reduce the risk of electric shock, burns, fire, etc.:

Always unplug this product from the electrical outlet immediately after use and before cleaning.

Please do not disassemble the product by yourself.

## WARNING

### Danger - To reduce risk of electric shock:

Do not operate under a blanket or a pillow; allow for normal heat dissipation so the product does not overheat, which could lead to fire, electric shock, or injury.

Do not unplug by pulling the power cord.

Only use this product with the power supply unit provided, and only with the voltage specified on the label.

Do not operate in places where corrosive products are used or where oxygen is administered.

Do not use in environments with high humidity (e.g.: bathroom) or high temperature, which could lead to electric shock or other safety hazards.

To avoid high temperatures, do not leave product in the car.

Do not use outdoors.

### Attention - Read all instructions carefully before using the product:

This product is a personal massage product, not a medical instrument, and it cannot be used for commercial activities.

Only use this product for the intended purpose as described in this manual.

Please switch the power off when the product is not in use.

The product should never be left unattended when plugged in. Unplug from outlet when not in use.

The product is equipped with a heating function. People with heat sensitivity should be careful when using.

Keep the cord away from heated surfaces, so as not to cause serious burns.

Please pay special attention to high wear components. Replace defective components immediately and / or keep the product out of use until repair.

Do not use the product if it has a damaged cord or plug, if it is not working normally, if it has been dropped or damaged, or if it has dropped into water.

Stop using the product immediately if you feel discomfort or pain during use.

Stop using the product if it is not working normally.

Do not use the product while driving, and observe all road regulations.

This product is not suitable for children, and children should not play with the product.

### If you have any of the following conditions, please consult doctor before using this product:

People with physical disabilities.

People with psychiatric issues or lack of experience and knowledge.

People with chronic illnesses.

People who are pregnant.

People who are menstruating.


People with physiological or other health issues.


People with pacemakers or concerns about their health.


## MASSAGE AREAS

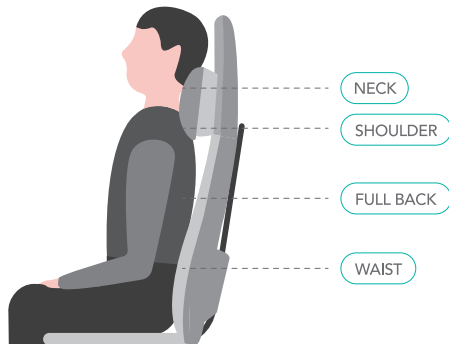
## PRODUCT FEATURES
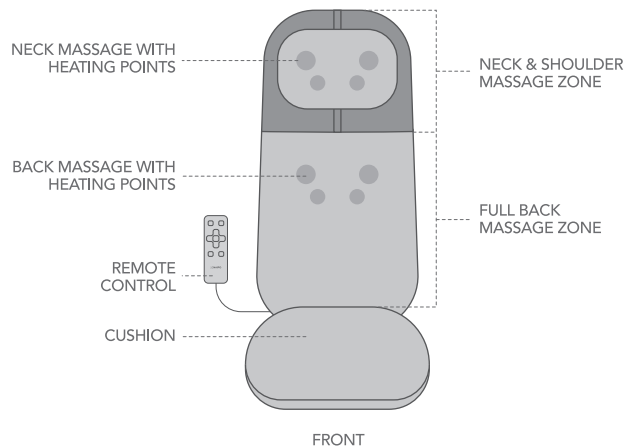
## IN THE BOX

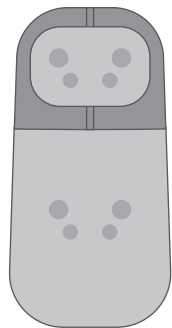
Massager

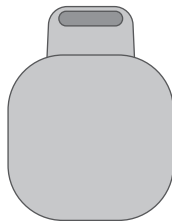Cushion

Wall Adapter
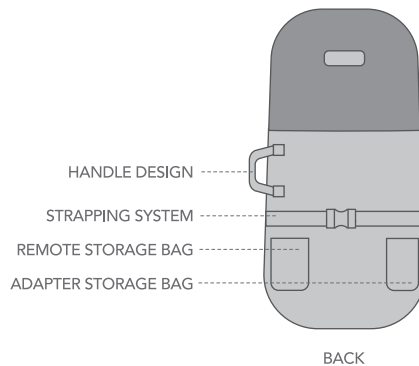
# OPERATION GUIDE

## How to use

Connect the product to power supply. Product will enter standby mode.

**Power:** While the product is in the standby mode, press the Button. The power indicator light will turn red, indicating default massage mode is on: Neck and Back Massage Level = "Low", Heat Level= "High".

**Timer:** While product is powered on, the default timing is 10min. Press the Button to select the massage duration: 10min / 20min / 30min. The product will shut down after the selected time has elapsed.

**Neck:** While the product is powered on, press the Button to start the neck pillow adjustment function: Long press the Up / Down to adjust the neck pillow up or down continuously. Release the button and the neck pillow will stop travelling.

**Back:** While the product is powered on, press the Button to start the back massage adjustment function or the fixed-point back massage mode: Long press the Up / Down to adjust the back massage up or down continuously. Release the button and the massage nodes will stop travelling, and proceed the fixed-point back massage.


**Auto:** While the product is powered on, the neck and back massage start automatically in the default massage mode. When using other massage modes, press the Button to switch back to default mode.

**Heat:** While the product is powered on, "High" Level will be turned on in the default mode. Press the Button to change the heat level. The cycle switches from High to Low to Off.

**Speed:** While the product is powered on, press the Button to change the speed of the back massage. The cycle switches from Low to Medium to High.
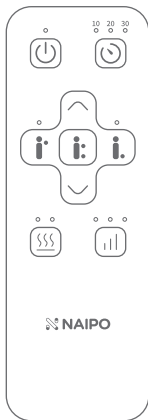
## Neck or Back Massage Function:

**Neck Massage only:** Long press for 2 seconds to turn off the neck massage function. This will not affect the back massage function. Press again to restart the neck massage function.

**Back Massage only:** Long press for 2 seconds to turn off the back massage function. This will not affect the neck massage function. Press again to restart the back massage function.

When the neck and back massage are operating separate massage functions, press to reset the neck and back massage functions.

## Heat Function:

Long press for 2 seconds to turn on the heat mode. The neck and back massage functions will stop. Press to toggle between "High" and "Low".

## In Heat Mode:

Press to switch back to the default massage mode.

Press to start the neck massage function, and press the Up / Down to adjust the position of the neck massage.

Press to start the fixed-point back massage, and press the Up / Down to adjust the position of the back massage.

## Please note during use:

- Please use it 30 minutes or less per day.
- Do not use while sleeping.

## CLEANING AND STORAGE

To prevent damage and discoloration, do not clean the product with corrosive detergent. Do not pierce or scratch the product with sharp objects.

Before cleaning, make sure that the product has been turned off and that the adapter has been pulled out from the outlet. Clean with a dry cloth.

Keep the product clean, and store it in a safe place. Keep away from environments with high humidity or temperature.

## MALFUNCTIONING AND MAINTENANCE

| PROBLEMS  | SOLUTIONS  |
|---|--|
| The product is not working properly                   | Connect the plug to the power supply, and check the wiring.  |
|   | The preset timer is finished. Press the power button again to restart.   |
|   | The product may have been in use for too long and overheated, activating the safety-oriented protection features. Turn off the power and wait about an hour for the temperature to return to normal. |
| Only heating without massage nodes working            |  Press <b>Auto Button</b> to start massage function.  |
| Neck massage can not reach the corresponding position |  Press <b>Neck Button</b> to adjust the neck massage pillow position.   |

**NOTE:** Do not disassemble or modify the product by yourself. NAIPO will not be responsible for any loss caused by unauthorized disassembly or modification of the product. This action may put the user in danger and / or render the product unusable.

Please contact NAIPO customer service for further assistance.

**NAIPO shall not be liable for costs, damages or repairs incurred as a result of, due to, or from:**

- The warranty period having expired.
- Unauthorized actions such as: product repair, misuse, abuse, modification, disassembly, etc.
- Artificial or accidental damage.
- Damaged internal and external packaging or missing accessories.
- Damage caused by force majeure.

## FCC COMPLIANCE

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions:

- (1) This device may not cause harmful interference;
- (2) This device must accept any interference received, including interference that may cause undesired operation.

The manufacturer is not responsible for any radio or TV interference caused by unauthorized modifications to this equipment. Such modifications could void the user authority to operate the equipment.

**NOTE:** This equipment has been tested and found to comply with the limits for Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and the receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

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